

Name:


Jagathpala Sneh

Date:

8/12/61

Experiment:

039

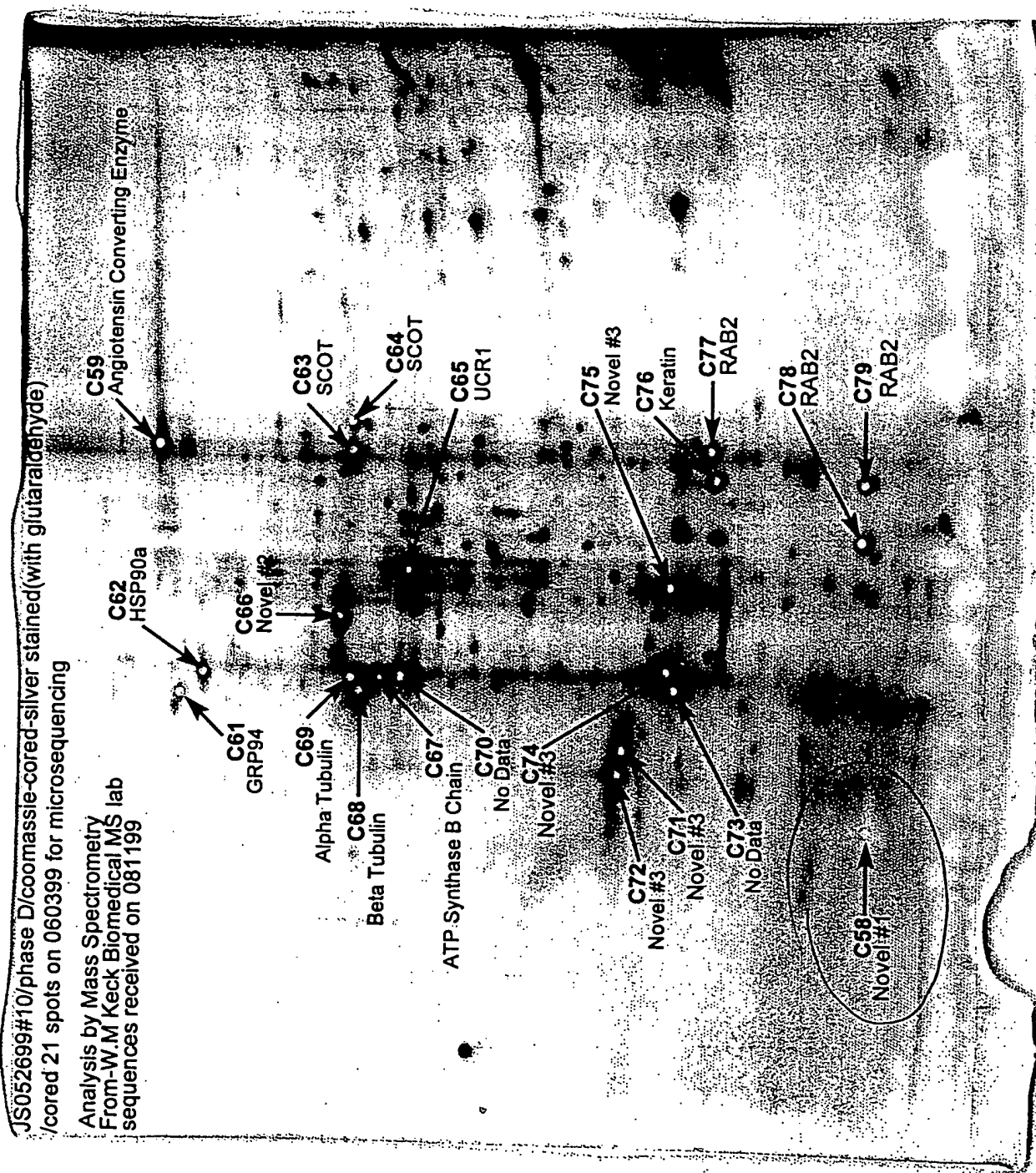

10/809,654

BEST AVAILABLE COPY

EXHIBIT 1

Name:

Jagathpala Srinivas

Date:

8/15/99

Experiment:

Report number: 400

Sequence Analysis of 22 2D Gel Bands.

8/11/99

Band C58. The peptides shown in Table 1 were detected in Band C58 (LB6-43-1). These peptides belong to Novel #1.

Table 1. Peptide sequences from Band C58 (LB6-43-1).

Peptide No.	Measured M W (M+H ⁺ , Da)	Peptide sequence by CAD ¹
1	1482.8 +2	ATSC ^a GLEEPVSYR
2	1499.4 +2	ATSC ^a (o)GLEEPVSYR
3	5033.8 +5	--- XSDSMEC ^a ---
4	5049.7 +5	--- XSDSM(o)EC ^a ---

GLEEPVSYR ~ 9mer

¹I and L cannot be distinguished by low energy CAD but are inferred by the database sequence, M(o) designates oxidized M, C is carbamidomethyl modified unless noted as C^a (acrylamide), _ designates a single unknown residue, - - - designates an unknown number of unknown residues.

Name: Jagathji la shah Date: 8/15/99
Experiment:

042

Nucleotide and deduced amino acid sequence of
Human Testis EST (AC # AA778671) which matched to
tryptic peptide obtained by Mass spectrometry of c

Soares Testis NHT Homo sapiens cDNA clone 1049023
mRNA sequence.
ACCESSION AA778671

```
1  GCACTGGTCCGGTCATCAACAAAGGCTGCCTGCGAGCCACCAGCTGCGGCCTTGAGGAAC
  T G P V I N K G C L R A T S C G L E E P 60
61  CCGTCAGCTACAGGGGCGTCACCTACAGCCTCACCACCAACTGCTGCACCGGCCCGCCTGT
  V S Y R G V T Y S L T T N C C T G R L C 120
121  GTAACAGAGCCCCGAGCAGCCAGACAGTGGGGGCCACCACCAGCCTGGCACTGGGGCTGG
  N R A P S S Q T V G A T T S L A L G L G 180
181  GTATGCTGCTTCCTCCACGTTTGCTGTGACGAACAGGGAGGACAGGGCCTGGGACTGTTC
  M L L P P R L L * P T G R T G P G T V L 240
241  TCCCAGATCCGCCACTCCCCATGTCCCCATGTCTTCCCCCACTAAATGGCCAGAGAGGC
  P D P P L P M S P C P S P T K W P E R P 300
301  CCTGGACAACCTCTTGCGGCCCTGGCTTCATCCCTTCTAAGGCTGTCCACCAGGAGCCCC
  W T T S C G P G F I P S K A V H Q E P G 360
361  GTGCTAGGGGAAGCATCCCCAGGCCTGACTGAGCGGCAGGGGAGCACGGCCCGTGGGTTT
  A R G S I P R P D * A A G E H G P W V * 420
421  GATTGTATTACTCTGTTCCACTGGTTCTAAGACGCAGAGCTTCTCACATCTCAATCAGGA
  L Y Y S V P L V L R R R A S H I S I R M 480
481  TGCTTCTCTCCATTGGTAGCACTTTAGAGTCCATGAAATATGGTAAAAAATATATATATA
  L L S I G S T L E S M K Y G K K Y I Y I 540
541  TCATAATAAATGACAGCTGATGTTCAAAA
  I I N D S * C S K 569
```

Name: Jagathpali Sheth Date: 05/26/9,
 Experiment: PCR to generate C58 - partial cDNA


043

PCR both primers for C58-EST
 using both forward and reverse primer.

Bottom:		Top:
3.025	3.8 bf	4.15-
2	4 dNTP	
2	0.2g	
1.25	CSPF (C58-F-EST)	
1.25	CSPR (C58-R-EST)	
0.475	H ₂ O	7.95
	cDNA	2
	polyase	0.5

PCR programme

- ① 94 2:30
- ② 94 1:30
- ③ 68 1:30
Δ-1.5/cycle
- ④ 68 2:30
- ⑤ 60-70 2 (1x)
- ⑥ 94 1:30
- ⑦ 50 1:30
- ⑧ 68 2:00
- ⑨ 60-70 6 (27x)
- ⑩ 68 18:00
- ⑪ 40 10

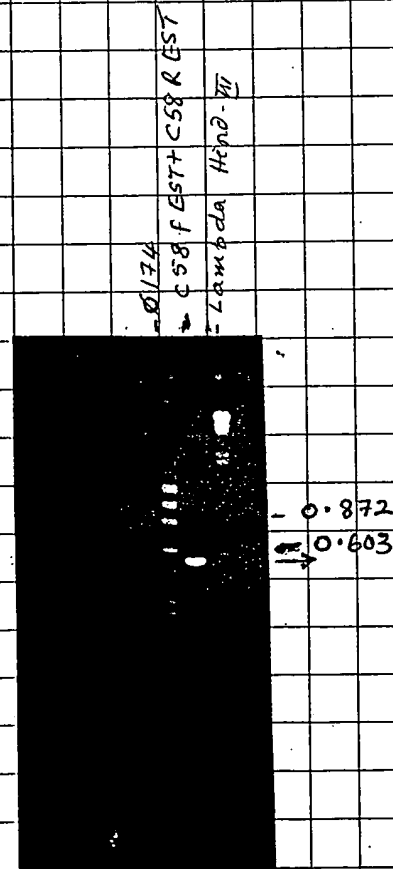

Result: Obtained a product around 530 bp. which matched to the expected product i.e. 519 bp

Name: Jagathpau Sheth Date: 9/7/99
Experiment:

048

The sequence for c58 est was
obtained from the sequencing lab.

Sequence of PCR-derived EST
9/7/99
partial sequence for c58

!!NA_SEQUENCE 1.0
Sequence of PCR-derived EST from 9/7/99
c58est.dna Length: 475 September 7, 1999 12:00 Type: N Check: 5379 ..

```
1 CTGCGGCTT GAGGAACCG TCAGCTACAG GGGCGTACC TACAGCCTCA
51 CCACCAACTG CTGCACCGGC CGCTGTGTA ACAGAGCCCC GAGCAGCCAG
101 ACAGTGGGGG CCACCACCAG CCTGGCACTG GGGCTGGGTA TGCTGCTTCC
151 TCCACGTTTG CTGTGACCAA CAGGGAGGAC AGGGCCTGGG ACTGTTCTCC
201 CAGATCCGCC ACTCCCCATG TCCCATGTC CTTCCCCAC TAAATGGCCA
251 GAGAGGCCCT GGACAACCTC TTGCGGCCCT GGCTTCATCC CTTCTAAGGC
301 TGTCCACCAG GAGCCCGGTG CTAGGGGAAG CATCCCCAGG CCTGACTGAG
351 CGGCAGGGGA GCACGGCCCG TGGGTTTGAT TGTATTACTC TGTTCCACTG
401 GTTCTAAGAC GCAGAGCTTC TCACATCTCA ATCAGGATGC TTCTCTCCAT
451 TGGTAGCACT TTAGAGTCCA TGA
```

Name: Jagadipali Sheth Date: 9/7/19. ^{testic cDNA}
Experiment: cloning of CS8 (Screening of Library) 050

A culture of K 802 strain host ~~cell~~ was made.

medium used : NZCYM medium.

20ml of NZCYM + 20% of 20% maltose soln
(Actual conc is 0.2% in the medium)

K 802 cell host strain taken from -70°C

with a sterile loop taken out and placed inside the medium.

Kept at 37°C - Shaker.

Name:

Jagathpala Sheth

Date:

09/08/99.


Experiment:

Cloning of c58

051

Titration of the λ DR2 Library.
(λ DR2 - human Testis cDNA Library)

- (1) The NZCYN medium was thawed using microwave.
- (2) About 20 μ l each of the medium was plated and poured on 5 plates and the cap was kept open (in the Sterile hood).
- (3) Mean time Take the culture of K802 left at 37°C previous day and take.
- (4) Take 1ml of ϕ SM buffer (buffer for λ DR2 i.e. phage buffer) in a tube and 9 μ l each to 5 tubes.

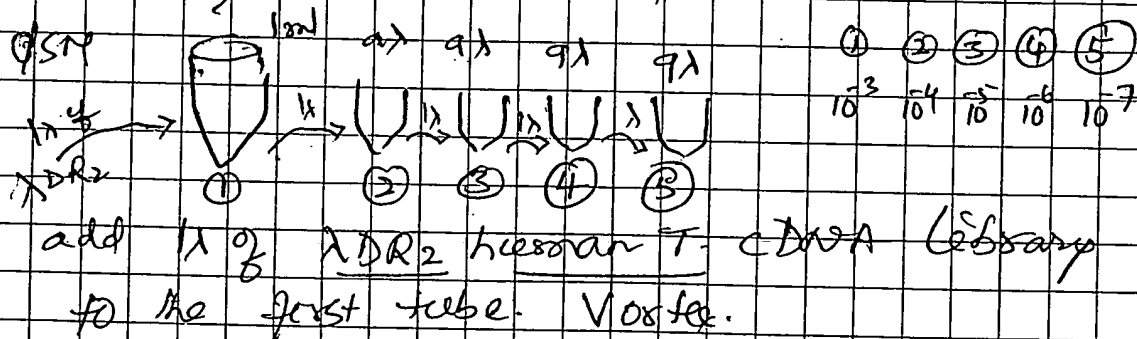

Take 1 μ l from tube 1 to tube 2, vortex and take 1 μ l from #2 & transfer to 3 and so on. vortex.

Take 1 μ l each from each tube and add to a 10 μ l tube (around bottom).

Name: Jagathpala Sheth

Date: 09/08/99

Experiment:


052

- (5) Add 75 μ l each of the K802 culture to all lobes. - wait for 20 minutes.
- (6) Melt the thaw NZCYM- agarose (0.7%) medium and allow it to come to $\approx 50^{\circ}\text{C}$ (for the top layer).
- (7) Keep a water bath at 37°C with a thermometer.
- (8) Keep the tubes at 37°C for 2 minutes.
- (9) Take $\approx 1/4$ out of the melted agar NZCYM agarose - in to the tubes ^{containing the} ~~pour~~ ^{pour} the contents from the tubes to the LB Agar plates, swirl the plates as you pour - Allow it to cool. for 10 min. to allow the inoculum to soak into agar.
- (10) Incubate plates at 37°C O/N.

Name: Sergathgalei Sheth Date: 09/09/19
 Experiment: Screening of Library

053

The plaques on the plates are counted.

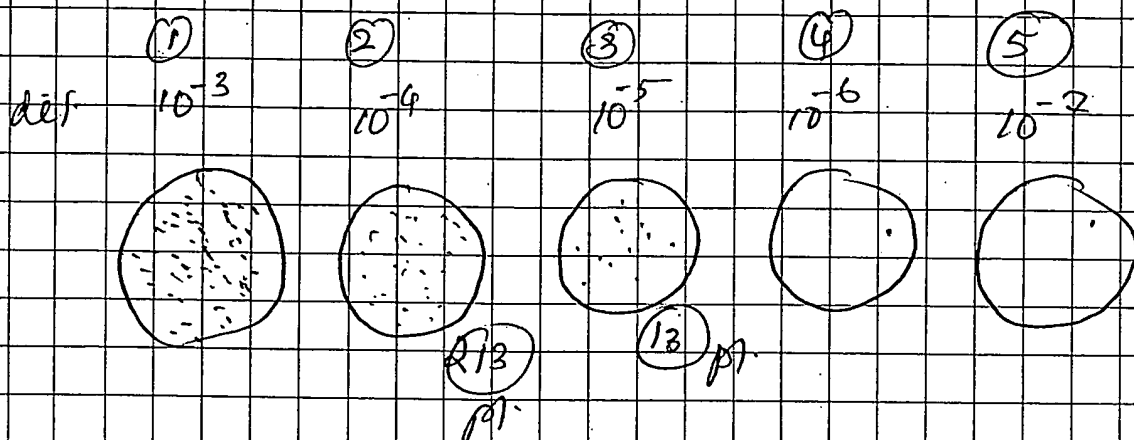

#(1) i.e. too many

#(2) 213×10^4 i.e. $2.13 \times 10^6 / \lambda$

#(3) 13×10^5 i.e. $1.3 \times 10^6 / \lambda$

average $\approx 1.7 \times 10^6 / \lambda$

average phage to be used for screening $\approx 40 \times 10^4$

phage
 λ dilute 100x

i.e. $\lambda \rightarrow 1.7 \times 10^7$

can take $\approx 2.5 \lambda$ i.e. ph

gives $\approx 50 \times 10^3$ phage

Name: Jagathpala Sheth

Date: 09/09/99

Experiment:

054

Infection of host strains

Poured 6 bigger plates with $\Delta 2$ cyp
poured 50 ml each (1.3% agar ^{medium})


Taken a small crystal of library ~~cell~~
from -70°C and the stock kept
back.

Take 1% \rightarrow det. 100 λ

205 λ \rightarrow should give \approx 50,000 phages

The bugs in 10 ml of $\Delta 2$ cyp with 2% maltose
— spun \rightarrow pellet taken and
resuspended in 10.5 ml of 10 mM
MgSO₄

Taken 600 λ each from 6 tubes
of bugs + 15 λ of phage.

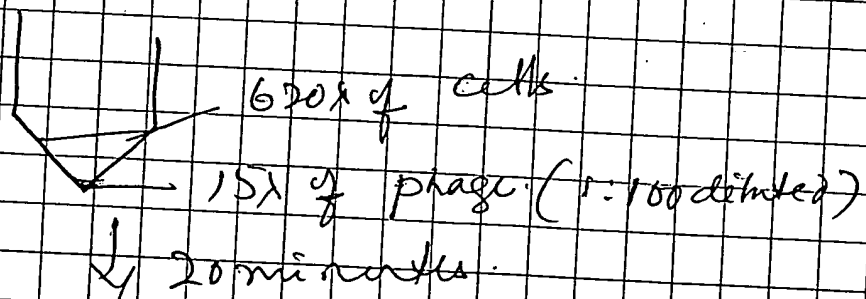

Take 100 λ each and
add to 6 tubes.
Clean line top agarose at 50-55 $^{\circ}\text{C}$.

10/809,654

EXHIBIT 10

Laboratory Research

National Brand

Name:

Jagathpala Sreth

Date:

09/09/99

Experiment:

Screening of Library

056

DNA Labelling

Protocol: Leinburg & Vogt Stein Method

To a sterile microfuge add:

C-58-EST

DNA

(= 50ng)

in 2 μ lH₂O3 μ l

50 min

oligo labelling

OLBf

10 μ l

bf

[α ³²P] dCTP5 μ l

Klenow

1.5 μ l

* After adding OLBf keep at -20 for a while.

Add 5 μ l of α ³²P dCTP and

1.5 μ l of Klenow. Incubate for a while and leave at 37°C.

Name:

Jyoti Pata Shethi

Date:

09/10/99.

Experiment:

cloning of c58 Contd. (Library Screening) 057

The plate - taken out from 37°C and chilled at 4°C .

membrane lifting.

① The nylon membranes - 6 of them numbered and 3 marks - were done at 3 corners - randomly.

② Membrane - placed on the plate carefully in one attempt. (Do not lift and change the position). - Leave for 2 min. (using tweezers)

③ Lift - make 5 marks with syringe needle. also make 1 mark on the side of the plate corresponding to position on filter. Take the membrane carefully and place it on a Whatman paper soaked

with denaturation soln - 5-10 min. Filter should be placed phase side up.

change positions in order ensure the complete immersion of the filter in the solution.

④ Place the membrane on a Whatman paper containing Neutralization buffer. change positions - ensure completely immersed and - 5-10 minutes



⑤ cross-linking:

- ① place inside on a Whatman
- ② press power on

Name:

Jagathisala Srethi

Date:

09/10/99

Experiment:

C58 - cloning - contd (screening of Library)

059

- * Open a food bag at one end.
- * Take the filter out using a folded whatman and put it to the bottom of the bag.
- * Seal ~~to~~ one side of the bag. - 2 seals.
- * Pour about 20 ml ^{formaldehyde} of the prehyb. solution
- * ^{same so ml for hybridization} Push the air bubbles out carefully.
- * Seal the top - 2 seals.
- * ~~Pour about 20 ml of the prehyb. soln.~~
- * Keep at 42°C - 3 hrs.

Purification of

Purification of the probe (DNA)

DNA purifying column - ~~end at the top~~

Remove the bottom cover. Cut the top off just below the matrix. Remove the plug off. ~~is~~ Take out the plunger. Insert into the

5 ml Syringe.

* ~~Load~~ equilibration of the column :- 5 ml of Elutip.

(low salt soln) = ~~slowly~~ Put the plunger and slowly bleed out the equilibration of. to a 15 ml tube.

* Take the labelled DNA (exude). Take \approx 900 ml of Elutip (low salt) - Elute out the equilibration buffer. Add one more ml of Elutip.

* Put the plunger and slowly push the plunger and get the unlabelled DNA to a 15 ml tube.

Add \approx 4 ml of the Elutip to the syringe. (Each time you reload the buffer disconnect the syringe, take the plunger out & then

Load sample)

* Disconnect the column. Connect to a

2 ml. fresh syringe (take the plunger out before connecting). Load \approx 1/2 ml of High salt solution.


Name: Jyothipala Shah Date: 09/10/99
Experiment: C58 Cloning - contd.

061

Replace the plunger and collect the
labeled DNA to a low microfuge
tube.

Hybridization.

- * Take out the membrane in the bag
Place a cut across the corner.
- * Pour off the soln to sink.
- * Take the purified probe and boil it for
5 min. (Open the tube in between (after
= 40 sec) and release press.)
- * Take 200 μ l of the hybridization buffer
(saved from earlier prehybridization
step) and add the labeled DNA
to it.
- * Pour this into the bag containing
membrane.
- * Carefully remove the air bubbles
out.
- * Seal safely - 2
- * Get the remaining bubbles to the
corner and seal again.
- * Clean all the areas.

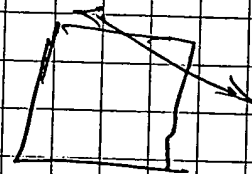

Name: Jagadishwala Sheth Date: 09/11/99

Experiment: C58- cloning- contd.

062

Washing of the membranes

- * Take out the bag, out the corner down to add container
- * Take out the membranes after cutting rough 3 sides
- * Place the membrane inside pour 200ml of the washing soln 1.

① Washing soln 1:

2X SSC made from 25X SSC in
SDS - 0.2% (200ml H₂O)

- Pour a small volume pour off after giving small short wash. Pour 200ml of solution and (solution at RT) and put the tray at 42°C - 20 min.
(It will come slowly to 42°C by 20 minutes)

② Washing Step 2:

0.2% SSC & 0.2% SDS - (200ml at 42°C)
prewarm the solution to 42°C.
incubate membranes - 20 min

③ Washing Step 3

200ml of 0.2X SSC & 0.02% SDS - 20 min
prewarm to 50°C (preferably 52°C)
incubate membranes - 20 min.

EXHIBIT 17

Laboratory Research

National CB and

Name: Jagathpali Sheth

Date:

09/13/99

Experiment: E58- cloning- contd.

063

Exposing the membranes

Take membranes in little 0.2 SSC and 0.2% SDS.

Take Cassette - mark

Place the int. screen on a flat surface on the bench place a s-wrap long enough. Place all the membranes in order.

fold the s-wrap. Place upside down. Fold the sides properly.

Place this on the cassette

the marked side up. (phase & side is down).

Take small piece of paper containing end pieces containing one or two dots and paste sandoship.

Place one Int. Screen on

the top.

Take to the client room place on a X-ray film & then another Int. screen - put at -70°C

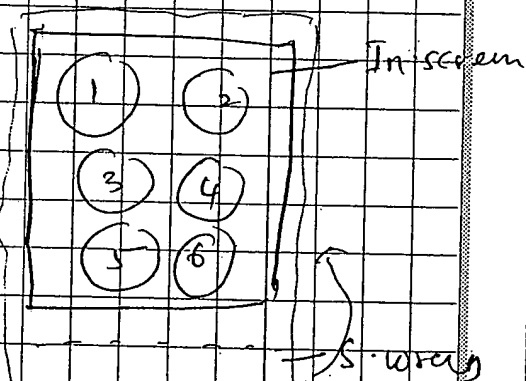


EXHIBIT 18

Name:

Jagadipali Sheth


Date:

9/1/99

Experiment:

064

9/1/99

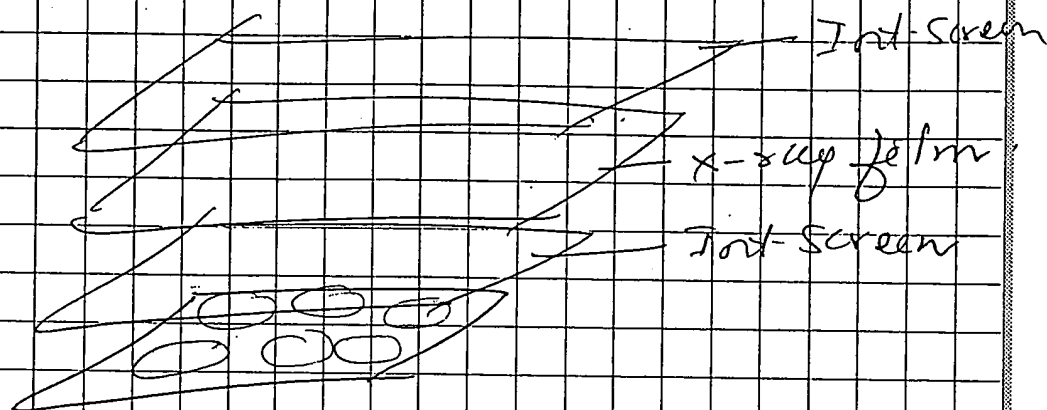

9/13/99

Exposed film taken out - One more film put in.


Align the film to the membrane and get to all the marks. Make an imprint on the X-ray film.

(preferably use diff. colours for different markings i.e. for periphery of the plates, side marks and 5 dot marks inside the membranes.

Mark the spots to be picked on the X-ray.

Name: Longathpala Shew Date: 09/13/99
Experiment: C58 cloning - contd

065

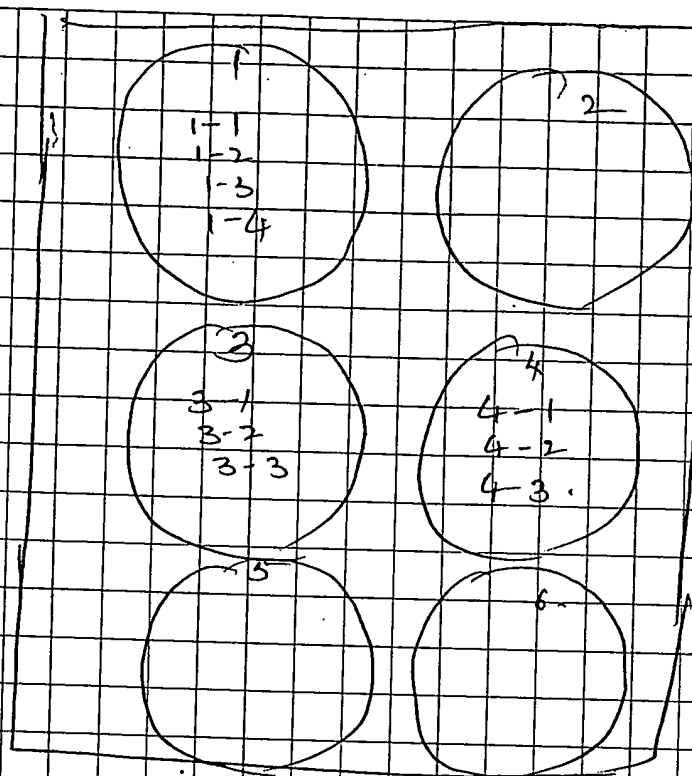

Decide about
the spots to be
picked.

~~Plus~~ Pipette 0.4 ml of DSM to 10 ml
tubes.

Aspirate the ~~sp~~ agar from the plate - shown
positive into tubes containing 0.4 ml
of SM.

put \approx 5% of Chloroform to each
tube. (increases in yield & also
sterilizes).

↓
put on a vortexing platform at
4°C for about 1-2 hrs

↓
Keep at 4°C till use.

Name: Jyothipala Sheth Date: 09-14-99.
Experiment: C58 cloning- contd.

066

Secondary Screening


Positive phage - taken out from 4°C

↓
Spin - 2 minutes

BSM → 90% to 10 tubes

1% of
+ phage sup. - 99%
PSN

↓
Vortex.

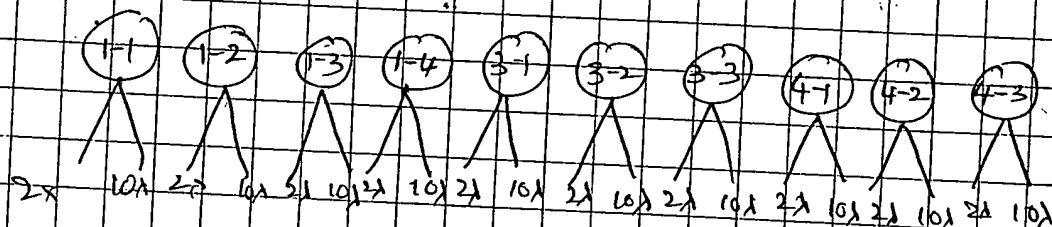

Mean time. 20 plates poured - NZCYM agar.
After solidifying, bottom - marked
with the ~~numbers~~ corresponding numbers.

NZCYM - agarose - poured - kept at 50°C .

Taken 2 tubes at a time containing
phage - kept at 50°C for 2 minutes.

↓
Add 4 ml of NZCYM - agarose

↓
poured a top layer on the plates
& allowed to solidify

↓
left at 37°C .

EXHIBIT 21

laboratory Research

National Brand

Name:

Jagathpala Suthi

Date:

09/14/79

Experiment:

C58 cloning control

067

DNA labelling:

50 μ g of C58 labelled as before.

8/

9/15/79

Secondary lifting

Plates taken out from 37°C.
In each pair the plate showing ≈ 200 phage selected

10 nylon (Bomall-8) - wrapped

A left was made as before

Denaturation (5-10 min)

NaOH (0.5M)
NaCl (1.0M)

Neutralization (5-10 min)

0.5M Tris
1.5M NaCl

Cross link

Drop the plates on membranes

Wash at 4°C with 2SSCA 0.2% SDS (30 min)

Name:

Jagathpala Shetty

Date:

09/05/99

Experiment:

C58 Cloning- contd.

068

C58
See

~~The~~ Prehybridization & Hybridization
Membranes

put in food bag (seal sides)

pour prehyb solution

3 hrs.

Purify the labelled DNA
using elutip. in 500s

500 μ g labelled DNA + 415 ml of
hyb solution

The bag opened & prehyb
soln. poured to sink

The hybridization done O/N
with the label + hyb soln.

EXHIBIT 23

Name:

Jagathpala Shetty

Date:

09/16/99

Experiment:

ESB - cloning - contd.

069

Washing of Membranes

① Discarded the lysis solution



I Wash 2x SSC, 0.2% SDS - 20 min. 30 → 42°C.
 (200 ml)

II Wash 0.2x SSC, 0.2% SDS - 20 min. 42°C.
 (200 ml)


III Wash 0.2x SSC, 0.2% SDS - 20 min. 42°C.
 (150 ml)



Membrane taken in 50 ml of 0.2x SSC & 0.2% SDS.



aligned on the Saran wrap.

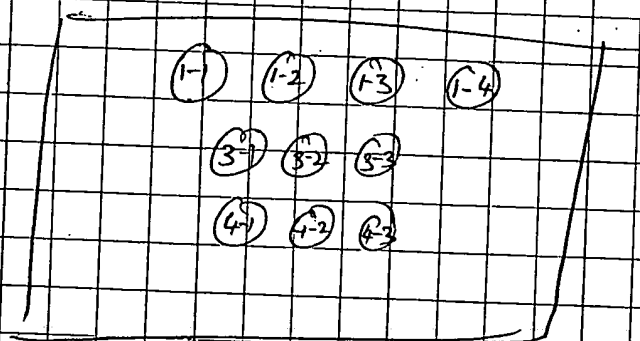

exposed at 11-45 AM.

09/17/99

Film developed and ~~to~~ marked respective to plates.

EXHIBIT 24

Name: J. Shetty

Date: 09-20-99

Experiment:


070

At 1, cells inoculated in 1 ml of LB + 1 μ g tetracycline
3 hrs. at 32°C shaker.

Spin the cells

Take pellet in 10 mM argSO₄ (1 ml)

Align the gel on the membrane and then
to the plate - mark the one true clone (isolated)
from the back of the plate.

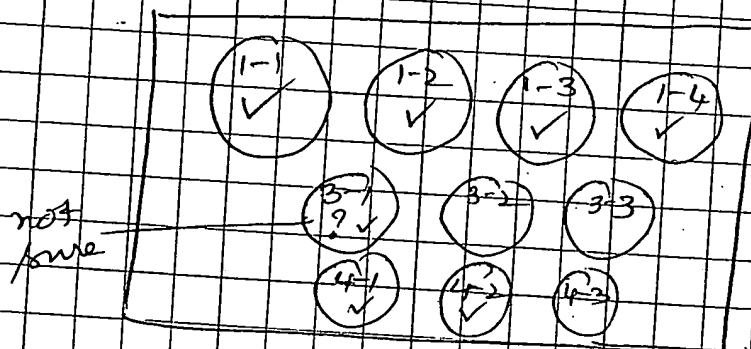

picked one
clone each
from the
marked ones.

into 1 ml tube
with 0.5 ml of DNA

5 \times of cell c₁

4°C shaker
1-2 hrs

Give a quick spin
to settle again

30 minutes } Take 15 \times from the sup

at 2 \times } 50 \times in 12.1.11

EXHIBIT 25

National Grand
atory Research

Name:

Jagathpala Shehi

Experiment:

c58- cloning- contd.

Date:

09-20-99

071

Continued from previous page

Add 50 λ of broth (LB) (recombination & circularization)
 1 hr. Shaker water bath at 32 $^{\circ}$ C

Add 2 μ l of 10 mM IPTG
 to induce replication of recombinant pDR₂

1 hr.
 Shaker water bath at 32 $^{\circ}$ C.

Add 1 λ of ~~500~~ 500 μ g/ml Carbenicillin
 & 1 λ of 1M Sod. citrate
 (for preferential existence of pDR₂ over pDR₁)
 32 $^{\circ}$ C for 1 hr.

5 λ \approx 50 λ
 Spread on LB-Agar plates.
 as follows.

20 λ of
 Carbenicillin
 500 μ g/ml

45 λ of Sod. citrate
 (1M)

5 λ or
 50 λ of
 Phage LAM1
 mix



Spread
 with Sterile
 Spreader

q/n: 37 $^{\circ}$ C

10/809,654

EXHIBIT 26

ratory Research

National Brand

FIR

CA

lin

Name:

Jagathpala Shethi

Experiment:

C58-cloning. Contd.

Date:

9/21/99

072

The plates observed and allowed to
grow to larger size at 37°C .

Left at RT for some time.

Inoculation of ~~to~~ to 3 ml LB cultures

Stock of 50 ml LB + 7510 ^{500 mg} Amp ^(50 mg/ml)
was made & divided 3 ml each / tube (15 ml tubes)

pick a single isolated colony using toothpick
choosing any one from a pair

inoculate to LB Amp.

Shake - water bath - 37°C o/n.
(for 3-1-2 colonies picked
i.e. 3-1a & 3-1b)

EXHIBIT 27

ory Research

linal® Brand

Name:

Jagadpala Sheth

Date:

09/22/99

Experiment:

Cloning of C58

073

O/N culture of A549 cells



Qiagen Kit isolation of DNA from plasmid

1. Cells pelleted at 2 steps.

Get - Take 1.5 ml into 1.5 tube - Spin (1 1/2 min)
discard supernatant, add another
1.5 ml and take the supernatant
using vacuum-dispenser.

Follow the Qiagen kit protocol to isolate DNA

① ~~dislodge pellet~~ Add 0.3 ml of Bf P1

dislodge pellet using P-200 pipettor

② Add 0.8 ml of Bf P2 - invert 4-6 ~~times~~ ^{times} - 5 minutes③ Add 0.3 ml of P3 - ~~gently~~ ^{gently} invert 4-6 ^{times} - keep ~~on~~ ^{at} ice - 5 minutes④ ~~spin~~ Spin - 10 minutes

⑤ Clean while set up the Qiagen column.

equilibrate the column with
1 ml of QBT

10/809,654

EXHIBIT 28

Name:

Jagathpala Sathi

Date:

09-22-99

Experiment:

Cloning of C58

074

- ⑤ Take the Supernatant Carefully from step ③ leaving the upper layer and the lower viscous pellet and load to the column Carefully.
- ⑥ Wash the column with 10ml x 4 (times) of solution QC. wait for last drop. Put the tube at the bottom of column.
- ⑦ Elute DNA with 0.8 ml of QF. wait till last drop.
- ⑧ Discard the column.
- ⑧ Add 0.56 ml of isopropanol.
- ⑨ Spin for 30 minutes, 12,000 rpm.
- ⑩ Take sup. with fine tipped pasteur pipette with a bulb.
(Make one fine tipped pasteur pipette)
- ⑪ ~~Give~~ Carefully load 200 μ l of chilled 70% ethanol and once again take the sup. off.
(* DO not disturb the pellet)
↓
Air dry

Important: -

Name:

Jegathpala Sneh

Date:

9-22-99

Experiment:

Cloning of c58

075

Dissolve DNA in 20 μ l of sterile water
 keep on shaker at 4 $^{\circ}$ C - 15 min

↓
 Shake again at \approx 20 $^{\circ}$ C with
 vortex mixer - 3-5 minutes

↓
 Give a quick spin.

Digestion of plasmid with BamHI
 and XbaI

BamHI

XbaI

isobut

(usually buffer conditions
 are different for 2 enzymes)

BamHI (Boehringer)

XbaI (Boehringer)

In this case
 same buffer used.

in 500 μ l tube

Vortex - give a quick spin Add last mix thoroughly & quick spin	{	DNA	- 2.5 μ l
		BSA	0.5 μ l
		10 \times BamHI B1	0.5 μ l
		XbaI	1.0 μ l
		BamHI	0.5 μ l

Name:

J. Shetty

Date:

9/23/99

Experiment:

Restriction digestion of DNA (plasmid)

077

Digestion of DNA - Sequential digestion

Cocktail for Xba I

9 λ of 10x Bf
 9 λ of 1mg/ml BSA
 19.5 λ of H₂O
 9 λ of Xba I

Add these,
 cool then
 add
 enzyme.

prepared
 for 18
 reactions.

Taken 1 λ of DNA + 4 λ of
 cocktail.

mixed with pipette tip.

37°C - 4-5 hrs

Bam-HI

Cocktail for Bam-HI

1.8 λ 5M NaCl - to bring the acetone conc
 to ~200 mM
 9 λ 10x Bam HI Bf
 9 λ Bam HI
 9 λ 1mg/ml BSA
 61.2 λ H₂O

prepared
 for 18 reactions

Added 5 λ each to
 tubes

= 37°C O/N

Name: J. Sheth

Experiment: Agarose gel electrophoresis of digested DNA

Date: 9/24/99

078

1.2% Agarose Gels

Lanes:

- ① 1-1 $\rightarrow M$
- ② 1-2
- ③ 1-3 - showed around 1 kb DNA band
- ④ 1-4
- ⑤ 3-1 $\rightarrow M$
- ⑥ 3-2 - showed around 900 b pair product
- ⑦ 4-1
- ⑧ 4-2 - showed around 1 kb product

①-② and ④-⑤

Given for Sequencing


DNA - 31

Forward primer 15A-23mer

Reverse primer 2A-21mer

150 11.5A
16A

1.078
872
603

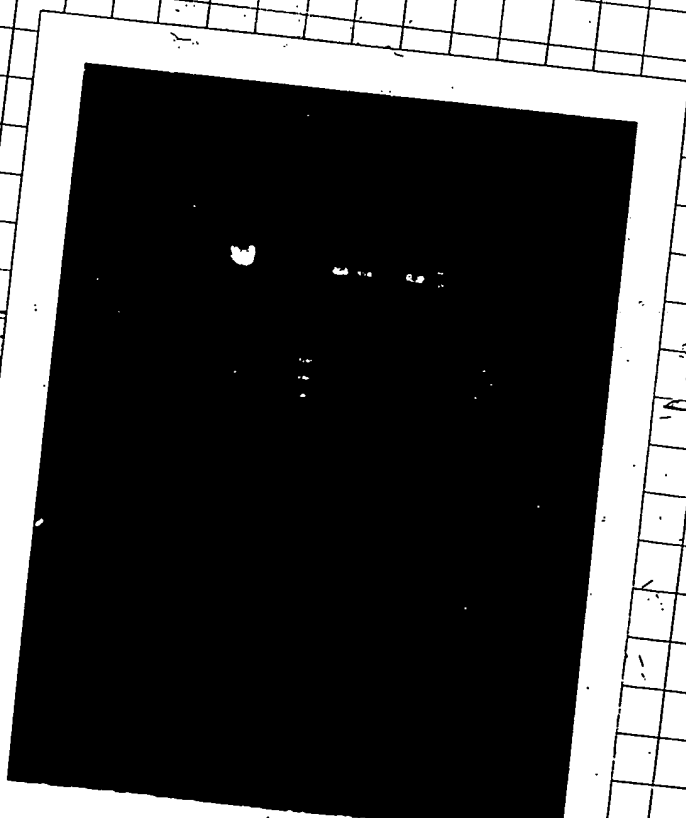

Research EXHIBIT 32

Research

NCAA

Cut the Sequence back.

Sequence - bad - as -

- Decided to give more DNA
for 1-2

DNA :	11.5 λ
Fer. pos m :	1.5 λ
H ₂ O	3.0 λ
	<u>16.1</u>

A culture of bacterial cells - with
clone - (1-2) and (4-2) (saved
earlier) - inoculated to LB
2% and culture with amp. and
Soc. Citrate. (10 mM) (75 μ g/ml)
O/A

Name:

J. Shetty

Date:

9/25/99

080

Experiment:

Plasmid isolation

Kidney - preparation of plasmid DNA

25 ml culture



Spun into 2 15 ml tubes

↓ spin - 3000 rpm.

~~pellet~~ Supernatant discard
completely

pellet



Processed for DNA isolation
using Qiagen Kit



pellet obtained at the
final step - ~~carefully~~
~~wash~~ - ~~add~~ carefully washed
with 2 ml of EtOH - ~~not~~ chilled

pellet dried completely



resuspended in 80 μ l of d.H₂O




saved in -70°C .

082

Sequence for 1-2 F - obtained -

Page
Tue, Oct 5, 1999 9
Mon, Oct 4, 1999 9
Spacing: 8.9

C T G G T A C C T A C A T G T G G C G A T G A C C G A G G A C T G C T T C A C A G G C C A C G G G G T C G C C C G G G A C T G T G T C C G G T C A T C A A C A A A G G C T G C T G C A G C C C A C C A G T G G G C C T T



 CACCTACAGCCTACCCACCACTGCTCCACCGGCGCTGTGTACAAAGCCCCGACAGCCAGACAGTGGGGGCACCACCACTCTGGCACTGGGGCTGGGATGCTGCTTCTCTCAG

0 590 600 610 620 630 640 650 660 670 680 690 700

EXHIBIT 35

Name: Jagadhpala Shetty Date: 10/6/91

Experiment: Sequence for C58

083

Sequence for 1-2 R & 4-2 R obtained.

However Sequences were bad.

They were resubmitted with a

request for $p(dt) > p(dt)_{20N}$

primer

However the sequence results of

Clone 1-2 F yielded enough

(good) Sequence to deduce the complete

Open reading frame for C58!

Name: Jagathpala S. 'h' Date: 10/6/99

084

Experiment:

Nucleotide and deduced
amino acid sequence for C58

Complete ORF of C58 contained 372 base pairs encoding 124 amino acids with a predicted Mol Wt. of 13 and a predicted pI of 5.5. Sequences of one of the tryptic peptides originating from the cored 2-D spot was found embedded in the ORF (Blue boxes).

GTCCCGGATCCGCGAGGGACGCAGGGCGTTGGGAACAGAGGACACTCCAGGCGCTGACCC
V P D P R G T Q G V G N R G H S R R * P -
TGGGAGGCCAGGACCAGGGCCAAAGTCCCGTGGGCAAGAGGAGTCCTCAGAGGTCCTTCA
W E A R T R A K V P W A R G V L R G P S -
TTCAGCGGTTCCGGGAGGTCTGGGAAGCCACGGCCTGGCTGGGGCAGGGTCAACGCCGC
F S G S G R S G K P T A W L G Q G Q R R -
CAGGCCGCCATGGTCCTGTGCTGGCTGCTGCTTCTGGTGATGGCTCTGCCCCAGGCACG
Q A A M ¹ V L C W L L L L V M A L P P G T -
ACGGGCGTCAAGGACTGCGTCTTCTGTGAGCTACCGACTCCATGCAGTGTCTTGGTACC
T G V K D C V F C E L T D S M Q C P G T -
TACATGCACTGTGGCGATGACGAGGACTGCTTCACAGGCCACGGGGTCGCCCCGGGCACT
Y M H C G D D E D C F T G H G V A P G T -
GGTCCGGTCATCAACAAAGGCTGCCTGCGAGCCACCAGCTGCGGCCTTGAGGAACCCGTC
G P V I N K G C L R A T S C G L E E P V -
AGCTACAGGGGCGTCACCTACAGCCTCACCAACCACTGCTGCACCGGCCGCCTGTGTAAC
S Y R G V T Y S L T T N C C T G R L C N -
AGAGCCCCGAGCAGCCAGACAGTGGGGGCCACCACCAGCCTGGCACTGGGGCTGGGTATG
R A P S S Q T V G A T T S L A L G L G M -
CTGCTTCCTCCACGTTTGCTGTGACCAACAGGGAGGACAGGGCCTGGGACTGTTCTTCCA
L L P P R L L ¹²⁴ P T G R T G P G T V L P -

Important: Place card under blue copy.

EXHIBIT 37

Name:

Jagathpala Shetty

Date: 2/10/99

Experiment:

Recombinant expression of C58

085

Primers ordered for the generation
of C58 - O55 - DNA - with Xho and
Nco site on either side to be
ligated to a PET 20 vector

Name:

Jagathpali Shethi

Date: 11/2/99

Experiment:

PCR to generate C58-Complete ORF

086

PCR reaction

use C58 PET primers.

Bottom~~Top~~Top

3.025

3.325

4.55

2

4 dNTP

2

Mg

1.25 (pm/A)

GSP-F' (C58 PET F)

1.25 (20 pm/A)

GSP-R' (C58 PET R)

0.425

1.20

7.95

1 DNA

2

poly dle

0.5

① C58 PET-R-60 pm/A

② C58 PET-R-20 pm/A

PCR programme (JBL)

① 94°C 2:00

② 94°C 1:30

③ 72°C 2:30

Δ-1/cycle

11 times

④ Go to ②

370

⑤ 94°C 1:30

⑥ 60°C 1:30

⑦ 72°C 2:00

27x

⑧ Go to ⑤

⑨ 72°C 18:00

⑩ 4°C ∞

⑪ END

Result: Gave the expected size product



Important: Place card under blue copy.

EXHIBIT 39

Name:

Jagathpala Sheth

Date:

11/16/99

Experiment:

089

Digestion of C58-PET-DNA with $XhoI$ and $NcoI$ endonucleases.

DNA received in 90% \rightarrow auto-evaporated to 15%.

Digested with $XhoI$ & $NcoI$ as follows:

(omega) bf	DNA	15%
	D (10%)	20.5%
	XH₂O	3.5%
Bovinger	$XhoI$	2%
NEB.	$NcoI$	2%
		25%

$\rightarrow 37^{\circ}C$ O/N
with DNA.

11/17

Amplification of DNA by gel electrophoresis

loads all 25% + 4% of loading bf.

used 5 wells (covered with tape) 5 wells.

DNA preserved in $\approx 80\%$ of d-H₂O

↓
De salted using Ambion x 2 times

↓
recovered in 60%

↓
Amplified.

Important:

EXHIBIT 40

Laboratory Research

National Brand

Name:

Jagadpala Shetty

Date:

11/17/99

Experiment:

090

Samples: ① 58 μ g DNA + 1 μ g loading bf.


② 31 μ g DNA + 2 μ g

③ 18 μ g DNA + 4 μ g

④ 0.5 μ g marker

⑤ 0.1 μ g marker

⑥ 1.5 μ g marker

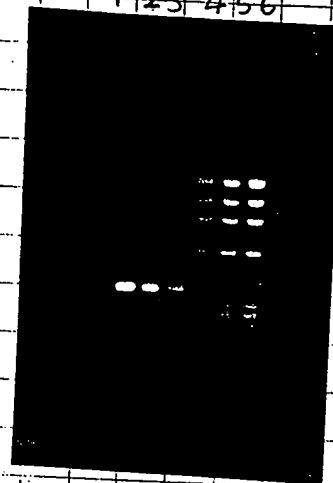

Actual amount of DNA: $\frac{125 \times 1 \times 603}{5.386} \times 1$

$= 13.99 \text{ ng}/\lambda$

Total vol: 50 μ l. i.e.:

699.73 ng

Name:

Jagatpala Sheth

Date: 11/18/99

Experiment:

091

Ligation

NCO/ Xho cut {

PET C58 -	3 λ
PET 286+	2 λ
10x lig-bf	2 λ (also contain ATP)
H ₂ O	12.5 λ
Ligase	0.5 λ
	20 λ

↓ 14°C
 O/N

After joining

DNA, λ & vector
 warmed to 50°C in
 water bath and for
 30 sec. and cooled
 at 25°C.

Then added lig-bf
 mix thoroughly and
 finally add ligase on
 ice.

Culture of host strain bacteria - Novo Blue.
~~DE3BL-21~~

1ml of LB + spec of strain - O/N 37°C.

11/19/99

Preparation of Competent cells
and Transformation of DNA to
host strains

① Culture diluted 5 times and checked
 O.D. (622nm) = 0.8
 Novabine = 0.8

② Diluted the culture back down to

10/809,654

EXHIBIT 42

Laboratory Research



Name:

Tegathala Shetty

Date:

11/11/97

Experiment:

092

0.1 OD in 1.25 ml LB + 12.5 ml MgCl₂/SO₄

ie: 170x ~~1.25~~ 25 g culture used.

Grown to 0.55 OD at 37°C Shaking

Centrifuged remove supernatant

Redissolved in 0.4 ml TFB (from NJW)
and keep in ice - 10'

Centrifuged, dissolve - 100 μl of TFB

Add 3.5 ml DMSO (from NJW)

Keep in ice 10'

Add 3.5 μl DMSO again
Keep in ice 10'

Add 10x each of ligation mixture
& kept in ice - 30'

Given a heat shock @ 42°C for 90 sec.

Kept in ice 2'

Added 300 μl LB + MgCl₂/SO₄ + Glucose
3x (ie: 8x for 1 ml)

Shaken at 37°C - 1 hr.

Name: Jyothipala Sneh

Date: 11/19/99


Experiment:

093

Plating.

Plating was done on LB-agar plates.

from each tube 3 plates were plated for BL-1 & number each out 40x, 360x & 45x 100µl of canamycin as the selection

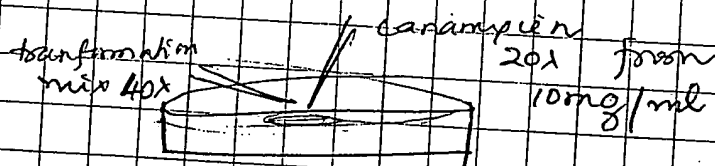

37°C O/N

11/20/99

One colony picked from plates 40x 360x ① & ② from each strain and a 3ml O/N culture made in LB + tetracycline + canamycin 10µg/ml

Isolation of Plasmid DNA

Isolation made by following the protocol in the Qiagen kit for miniprep.

- ① 3ml of culture - centrifuged in 1.5ml microfuge tubes at 2 steps.

- ↓
② Add 0.3ml of bf. P1 to the pellet dislodge the pellet with pipette.

- ↓
③ Add 0.3ml of bf P2 - invert 4-6 times at 4°C - Sit - 5 min

- ↓
④ Take P3 from 4°C and add 0.3ml to tube and invert 4-6 times and place it on ice - 5 minutes

- ↓ Spin 10 min
⑤ Meanwhile set up the Qiagen column. Equilibrate the column with QBT

↓
Take the supernatant from step 4 carefully and load to the column.

- ↓
⑥ Wash the column with 1ml x 4 of QC. wait till last drop drops off.

EXHIBIT 45

10/809,654

Name:

Jagathpala Shetty

Experiment:

Date:

11/21/99

095

⑦ Elute DNA in 0.08 ml of QF
wash till the last drop

↓

⑧ Add 0.50 ml of isopropanol

↓

⑨ Spin for 30 minutes at 10,000 RPM

↓

⑩ Remove sup. with Gene Lepper
piston pipette

⑪ Carefully wash the pellet with
200 μ l of chilled 70% ethanol.

↓

⑫ Air dry.

11/22/99

Resuspend the DNA in 20 μ l
each of sterile water.
mix at 4°C for 15-20 min

EXHIBIT 46

Name: Jagathpala 'sheli' Date: 11/22/99

Experiment:

096

Digestion of plasmid DNA with XhoI and NotI


	H ₂ O		
	DNA	2.5 μ	
(Boehringer)	XhoI	1 μ	} - 37°C O/N
(NBI)	NotI	1 μ	
(Pobonaga)	Bf D (10x)	0.5 μ	
		<u>5 μ</u>	

11/23/99

2% agarose gel electrophoresis of
digested DNA

- ① Marker ϕ
- ② B2-21-①
- ③ B2-21-②
- ④ Nov-B1-①
- ⑤ Nov-B1-②
- ⑥ Marker \times HindIII
- ⑦ Marker ϕ

① ② ③ ④ ⑤ ⑥ ⑦

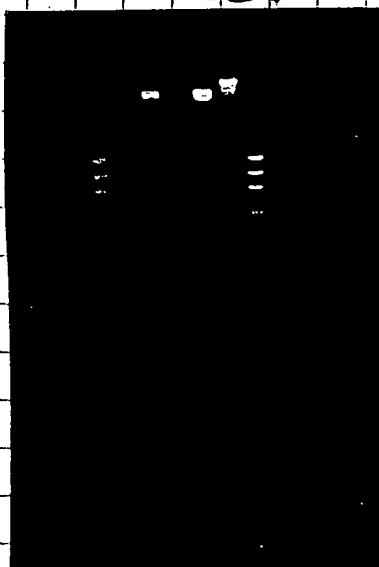

clone #④ (Nov-B1-②) gave the right size insert.
A Colony stock of the same - done

Name: Jagathpaal Shetty Date: 11-2-99
Experiment: Sequencing of the vector.

097

DATA from

clone #4.

Wara. Blue - ②

=

was

given

for

sequencing.

①

DATA :

~~8~~

8 λ

T7 terminator

2 λ

(5 pmoles/ λ)

H₂O

6 λ

16 λ

②

DATA :

8 λ

T7 promoter

—

requested from Bank

H₂O

4 λ

12 λ

098

With 2 enzymes: NCOI XHOI

November 29, 1999 14:35 ..

[illegible]

C58- is successfully ligated
to the pET28b vector

10/809,654

FYHIBIT 29

Name: Jagathpala Sathy Date: 11-25-97
Experiment: _____

099


A 0.1% culture from pET28b-C58-NvaBla (#4)
was made in tubes
(3 ml each)
↓
plasmid DNA isolated.

11/26/99

11/29/99 A 2% agarose gel run.


- ① Marker
- ② Tube # 1 from pET 28b-C58-NvaBla #4
- ③ Tube # 2 from " " "
- ④ Marker

1 2 3 4

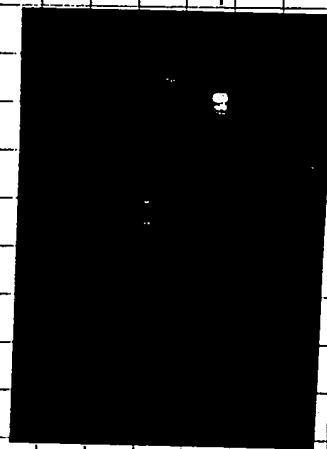


The host strains bearing the plasmid
-gene had 0.7 was diluted to
streak a plate, force a single colony
and make a glycerol stock of
the construct.

Fig. 8. Proposed Architecture of C58

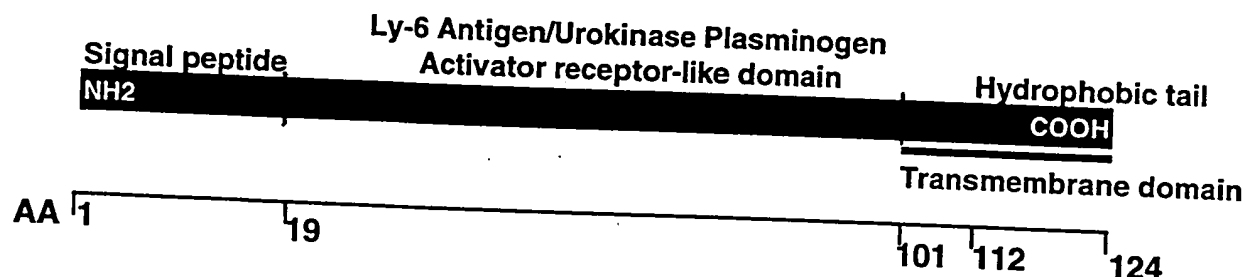

C58 is GPI anchored - It has a signal peptide, a transmembrane domain and a transamidase cleavage site!

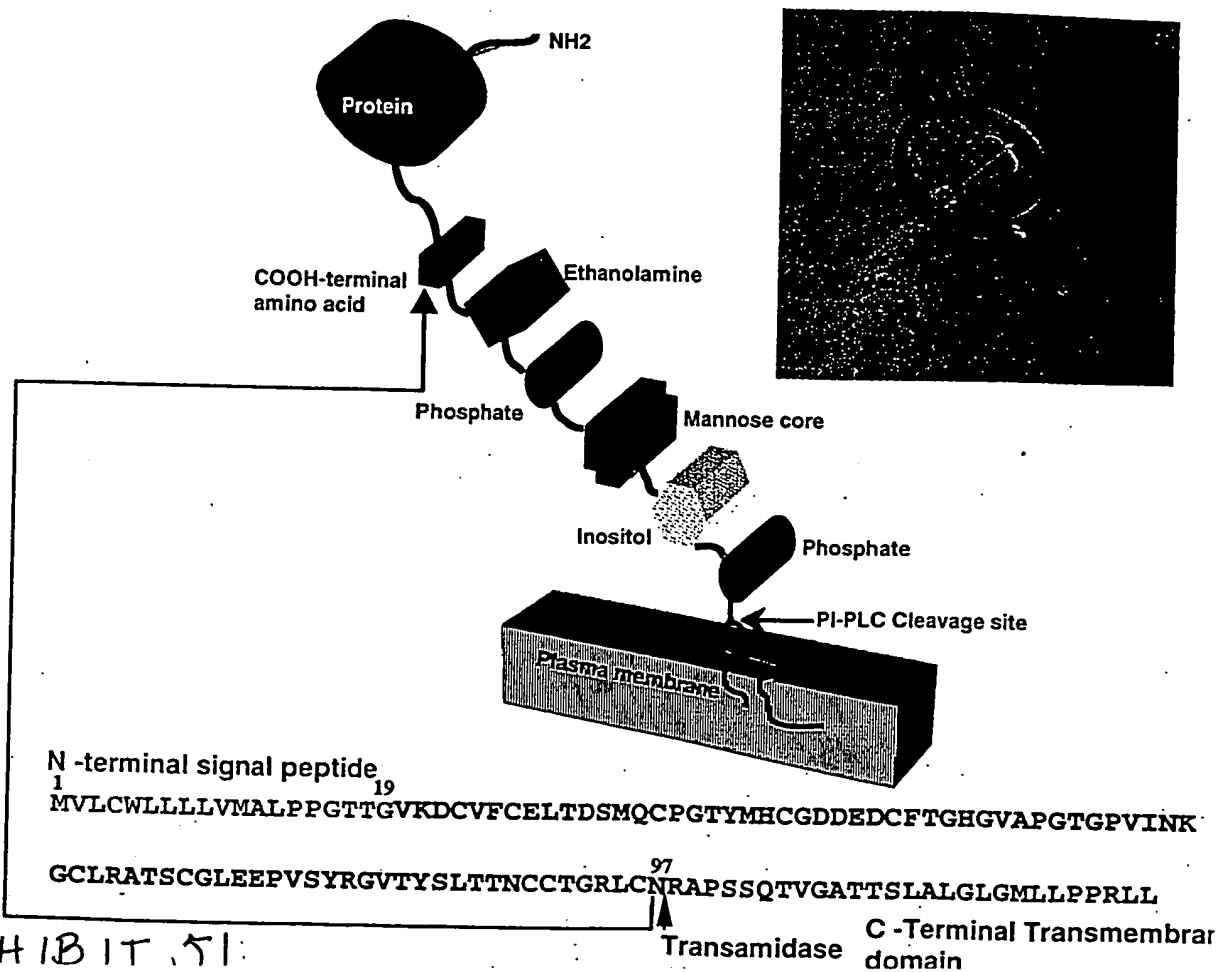

Name: Jagadevi - Shetty Date: 11/23/72
Experiment: Sequence analysis of C58

Sequence alignment of C58
with other Ly6/UPAR family
members.

C58 (24-98): VFCELTDSMQCPGTYMHCGDDEDCFTGEGVAPGTGPVIN---KGCLRATSCGLEEPVSIRGYTYSLTTCCTGRLCNRA
CD59-AOTTR (12-126): CPYPTTQ---CTMTTNTCTSNLDSCLIARA-GSRVYYR-----CWKFEDCTFSRYSNQLSEN-ELKYCCCKKNCNPN
CD59-CALSQ (12-126): CPYSTAR---CTTTTNTCTSNLDSCLIARA-GLRVYYR-----CWKFEDCTFRQLSNQLSEN-ELKYHCCCKKNCNPN
CD59-SAISC (12-128): CPLPTMESMECTASTNCTSNLDSCLIARA-GSGVYYR-----CWKFDDCSFKRISNQLSET-QLKYHCCCKKNCNPN
CD59-CERAE (12-126): CPNPTTD---CKTAINCSSGFDTCCLIARA-GLQVYNQ-----CWKFANCNFNNDISTLLKES-ELQYPCCKKDLNPN
CD59-PAPSP (12-124): CPNPTTN---CKTAINCSSGFDTCCLIARA-GLQVYNQ-----CWKFANCNFNNDISTLLKES-ELQYPCCKKDLNPN
CD59-HUMAN (12-126): CPNPTAD---CKTAVNCSSDFDACLIARA-GSGVYYR-----CWKFKCSFKRISNQLSET-QLKYHCCCKKNCNPN
CD59-HSVSA (7-117): CSHSTMQ---CTTSTCTSNLDSCLIARA-GSGVYYR-----CWKFDECNFDFISRLAEK-KLYNCCCKKDLNPN
CD59-PIG (12-123): CINPAGS---CTTAMNCSSHNQDACIFVEAVPPKTYQ-----CWRFSDCNAKPIILSRLEIA-NVQYRCCQADLCNPN
CD59-RAT: (9-120): CLDPV-SS---CKTNSTCSPNLDACLVAVS-GKQVYQ-----CWRFSDCNAKPIILSRLEIA-NVQYRCCQADLCNPN
LYGA-MOUSE (2-134): CYGVVFET-SCP-SITCPYPDGVCVTQEAIVIVDSQTRKVKNNLCPLICPPNIESMEILGTV-NVNTSCCKEDLCNA-
LYGP-MOUSE (11-107): CLGVSIGI-ACK-SITCPYPDAVCISQVELIVDSQRRKVKNNLCPLICPPNIESMEILGTV-NVNTSCCKEDLCNA-
LYGC-MOUSE (2-131): CYGVPIET-SCP-AVTCRASDGFCIAQNIELIEDSQRRLKTRQCLSPCPAGVP---IKDPNI-RERTSCCSEDLCNA-
LYGE-MOUSE (11-107): CTDQKNNI-NCLWPYSCQEKDHYCITLSAAAGFGN-YNLGYTLNKGCSPICPSENYNLNLGYA-SYNSYCCQSSFCNPN
E48A-HUMAN (21-93): CTSSSN---CKHSYYCPASSRFPCKTNTYEPLRGNLYK---KCAESCTPSYTLQGGYSSG-TSSTQCCQEDLCNPN
THYB-MOUSE (3-117): CTNSAN---CKNPQYCPNPFYFCKTNTYEPLRGNLYK---KCAESCTPSYTLQGGYSSG-TSSTQCCQEDLCNPN
UPAR-RAT (17-132): CESNQD---CLYECCALGQ---DLCRTTYLREWQDDRELEYTRGCAHSEKTNRTMSYRMGSMIISLTETCATNLCNPN
UPAR-MOUSE (14-131): CESNQD---CLYECCALGQ---DLCRTTYLREWQDDRELEYTRGCAHSEKTNRTMSYRMGSMIISLTETCATNLCNPN
UPAR-HUMAN (14-129): CKTNGD---CRYECCALGQ---DLCRTTYLREWQDDRELEYTRGCAHSEKTNRTMSYRMGSMIISLTETCATNLCNPN
UPAR-BOVIN (5-127): CENTTS---CSYEECTPGQ---DLCRTTYLSYWEGGNEMNYRKGCTHPDKTNRSMSTRAADQIITLSETYCRSDLCNPN

Name: Jagathjini Sheth

Date: 11/29/99

Experiment: Bacterial expression of C58

01

Bacterial cells (NOVA 022VE) containing
the construct pET 28b - C58 (#4)
was streaked on a agar plate (LB)

11/30

picked a single colony and inoculated
to 1ml LB broth

↓
a glycerol stock made
(1ml of culture + 150μl of 100% glycerol)

Protein Expression

0.5 of the culture from
above taken - inoculated
to 2ml LB culture medium
+
kanamycin - 10μg/ml

↓
gallon to ~ 0.5 OD
~~leave at~~ ↓
4°C O/N
↓

12/1/99

cultures from above inoculated
to 20ml culture (LB + Kanamycin)

MA
EX
Name: Jagadala Sheth
Experiment:

Date: 1/1/99

0

Expression - continued
Expt.

20 ml culture

checked
O.D.
600nm
(200 + 800 nm f)
LB

0.07 O.D.

induced
with 1mM IPTG

Stock 200mg/ml (840mM)

4 samples with 0.5 O.D/ml
saved, at 0 time
of induction)

O.D. after 2 hrs = 2.0
sample collected
after 2 hrs after
induction.
(0.5 O.D/ml
= 4 samples)

O.D. after 3 hrs = 2.8

4 0.5 O.D/ml samples
saved

Kept on ice.

Centrifuge

save pellet
at freezer

control

20 ml culture

0.07 O.D.

not induced

0.5 O.D/ml
samples (4)
saved

sample collected
after 2 hrs
after
0.5 O.D/ml
= 4 samples

O.D. after 3 hrs = 2.8

4 0.5 O.D
samples saved

Kept on ice

Centrifuge

save pellet at
freezer

10/809,654

Important: Place card under blue copy.

EXHIBIT 54

Bacterial lysate preparation and electrophoresis

- ① total cell preparation
- ② soluble fraction
- ③ insoluble fraction

Total cell: 0.5 OD pellet + 20% of 0.5 ml
10 mM Tris, pH 8.0
+ 20% of sample buffer
Boil 40 70°C - 2 min
Centrifuge
Load everything.

② Preparation of soluble & insoluble fractions
Used Rong Buster - novagen
0.5 OD/ml - pellet
+ 20% of Rong Buster
↓ vortex
Shaker 10 min
↓ Centrifuge
pellet ← Supernatant → soluble fraction
add Rong Buster 20%
mix vortex
Add 200 µg/ml lysozyme
incubate 5 min
Sample Bt
Load

insoluble fraction - continued.

↓
Add 6 vols of 1:10 lmgbuster


↓ vortex.

centrifuge

↓
pellet + 1:10 lmgbuster

↓ centrifuge

↓
pellet resuspended in
10mM Tris + sample buffer.

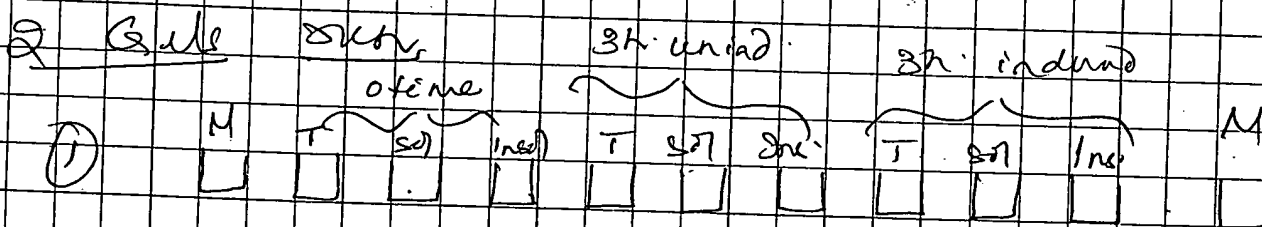

15% separating gel 4% stacking gel
run o/n at 150 mV

↓
Gel coomassie stained.

Name:

Tagahj. Ma Shells


Date:

11/14/99

Experiment:

05

Gel # 2.

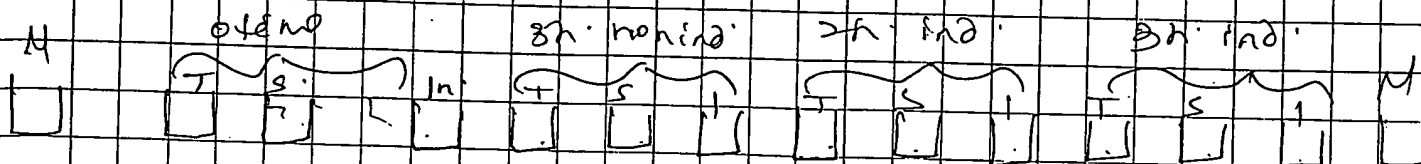

Gel. transferred to a nitrocellulose membrane

12/6/99

Western Blotting of the membrane

used 1:1000 dilution of Ni-NTA conjugate.


↓
developed by ECH
& TMB.

blot was prestained with
pierce's.
before probing.

Name: Jagathy la Sheth Date: 12/1/99
Experiment: _____

06

Gel # 1 : Coomassie stained
gel

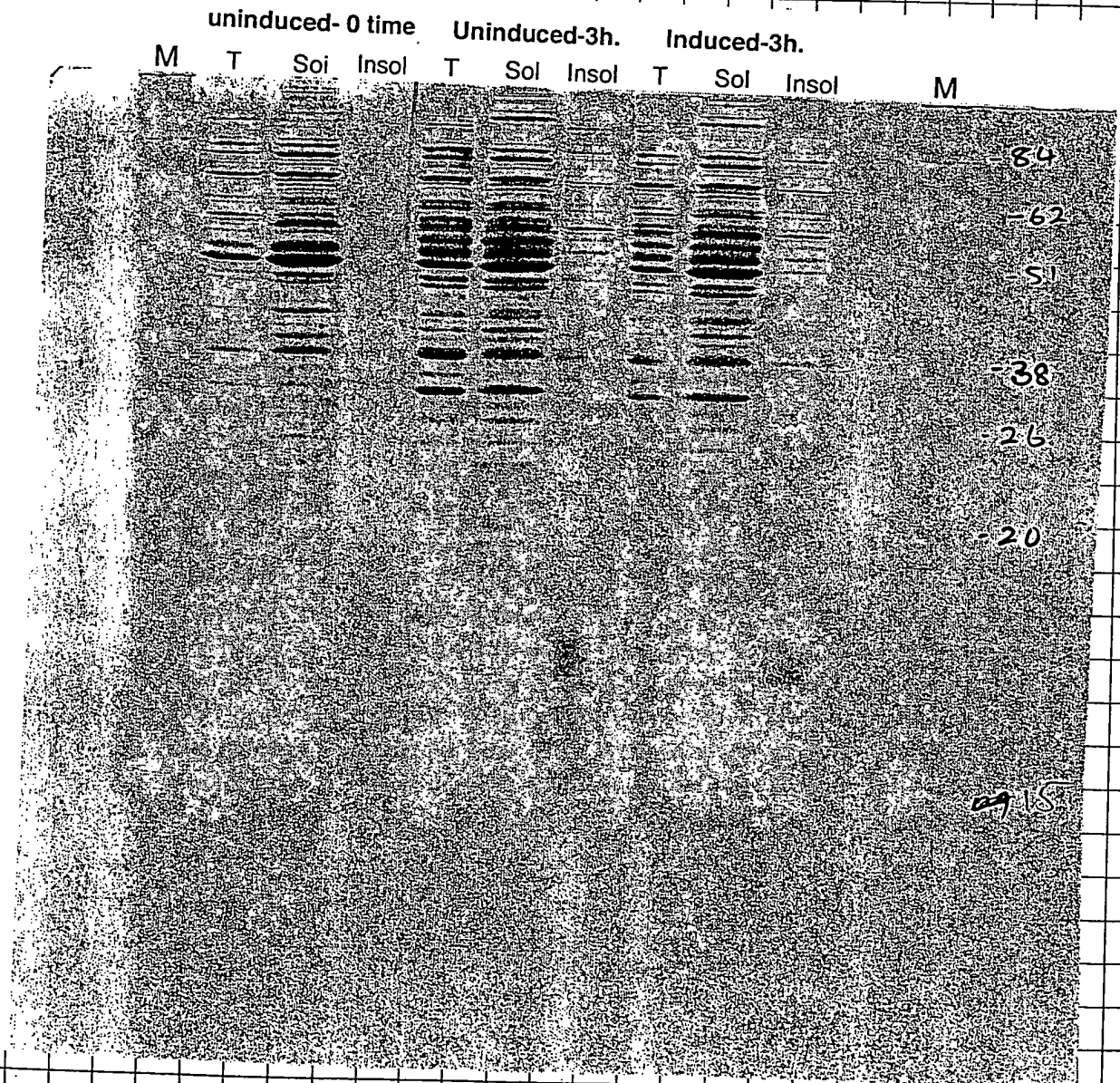

NO expression!

Important: Place card under blue copy. EXHIBIT 58.

Name: Jagathpaul Shetty

Date: 12/1/99

Experiment: C58 Rec. expression

09

- (1) Nova Blue - pet 28b - control - ~~used 7.6mg/l~~
 - (2) BL-21 - pet 28b - control - ~~used 7.6mg/l~~
 - (3) BL-21 - C-58-28b - transformed with C-58 plasmid DNA - 1x
- Plated 360x 2 40x

37°C o/n

12/11/99 Plates - examined and kept at 4°C

12/12/99

A single colony from one of the plates from each group - inoculated - 1ml o/n culture made



12/13/99 (1) Glycerol stock of all the 3 made.

(2) A 2ml culture for pet 28b-C58 & pet 28b - control made. till the OD reaches 1.0



kept at 4°C o/n

Name: Jagathpali Shetty

Date: 12/13/19

Experiment: C58- Rec. expression

10

documented 2 ml of inoculum from
control (empty vector) and C58+ vector-Novabne
to 20 ml culture.

At O.D = 0.7 added 10mM IPTG
to the culture.
Sample saved ~~at~~ before induction
(0.5 O.D/ml samples)

After 2 hrs. - samples saved
(0.5 O.D samples)

After 3 hrs. flasks taken out
chilled - ice.

0.5 O.D samples - aliquoted - centrifuged
Rest of the samples - centrifuged
and saved

↓
pellets saved at
-20°C

SDS-PAGE of the culture
Bacterial lysate

Gel: 15%

Sample preparation:

used bug buster - pellets dissolved
in 30 μ l of bug buster - vortexed - 5 min.
centrifuged \rightarrow sup. \rightarrow 30 μ l of sample at $\frac{90^\circ\text{C}}{2\text{ min}}$ - load
 \downarrow
pellet + 30 μ l bug buster

\downarrow
vortex
+ lysozyme 200 μ g/ml

\downarrow
incubate 5 min.

Add \downarrow 180 μ l of 1:10 bug buster

\downarrow
vortex

\downarrow Centrifuge 4°C 20 min

pellet \rightarrow add 200 μ l of
1:10 bug buster

\downarrow vortex
pellet \rightarrow add 200 μ l of
1:10 bug buster

\downarrow
pellet + add 200 μ l of
in this buffer 50 μ l

load

\uparrow
heat to 90°C
2 min

\uparrow
Add eq. volume
of sample \leftarrow in this buffer to

Name: Jagadipal Shetty

Date: 12/1/99

Experiment:

12

Controls: 200 ng of lysozyme in
50 mM Tris and sample buffer,
& sets of gels run.

Coomassie
Staining

Transferred to nitrocellulose
& probed with anti-NTA
(1:2000).

TMB

control
Unind Ind Unind Ind 2h Ind 3h Lysozyme (200ug)
M Insol Sol Insol Sol Insol Sol Insol Sol Insol Sol M

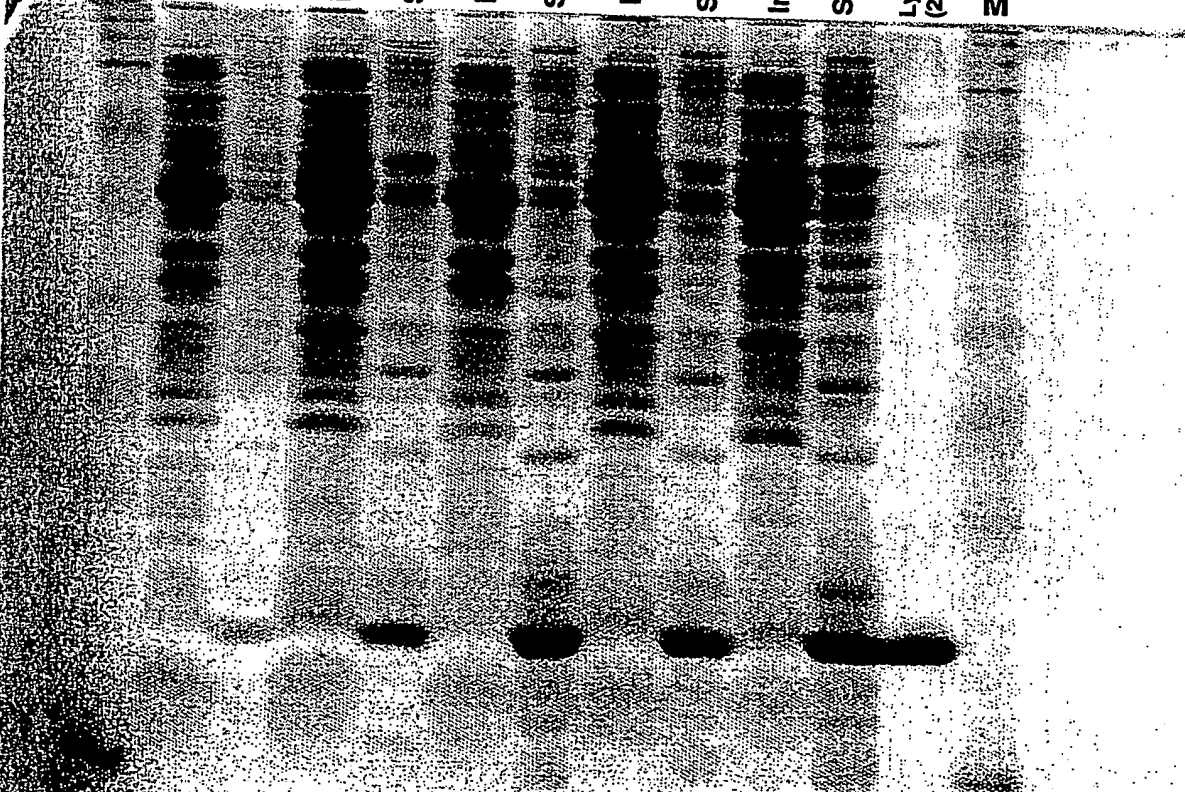


EXHIBIT 62

Name: Jagadhipali Sheth Date: 12/15/11
 Experiment: SAT- PAGE analysis

→ Marker (Gibco)
 — S } control (empty vector)
 — IN } unind.
 — S } control induced
 — IN }
 — S } c-58
 — IN } uninduced
 — S } c-58
 — IN } induced 2 hrs
 — 1 rest }
 — 50% } c-58
 — insoluble. induced 3 hrs
 — ~~rest~~ hypoxyme 200 µg
 — Marker

Name: Jagathpaul Sheth

Date: 12/15/99

Experiment: Northern Blot analysis of C58

14

Northern Blot analysis

Probe: c-58-ORF. The PCR product ~~is~~ of C58 ORF with Xho-I and Nco-I ^{BSA} on either side was cleaved ~~with~~ above enzyme and purified on agarose to clean off the end fragments.

Labeling of DNA (Vogelstein's method)

4x of DNA (13.9 ng/ μ l)
29.5x H₂O \rightarrow Boil - 5 min. in water bath
10x of DAB (from MJM)

\downarrow
Keep at -20°C for a while.

\downarrow
Add 5x of [α^{32} P] dCTP

Add 1.5x of Klenow polymerase

\downarrow
Incubate for a while

\downarrow
Keep at 37°C O/N

MTN® Blot User Manual

I. Introduction continued.

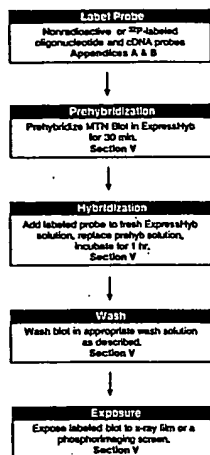


Figure 1. Overview of MTN Blot protocol. Use the β -actin probe to verify that hybridization procedures are working properly and to quantify results.

II. List of Components

Store unused MTN Blots at room temperature in a sealed plastic bag away from light.
 Store used MTN Blots at 4°C in a sealed plastic bag until needed.
 Store control probe at -20°C.

- 1 MTN Blot
- 100 ng Human β -actin cDNA control probe (2.0 kb) in 20 μ l of TE buffer (pH 7.5). Sufficient for 2-4 labeling experiments.
- 25 ml ExpressHyb™ Hybridization Solution

III. Additional Materials Required

- 20X SSC
 - 3 M NaCl
 - 0.3 M Sodium citrate (pH 7.0)
- Wash Solution 1
 - 2X SSC
 - 0.05% SDS
- Wash Solution 2
 - 0.1X SSC
 - 0.1% SDS
- Wash Solution 3
 - 2X SSC
 - 0.1% SDS

CLONTECH Laboratories, Inc.

www.clontech.com

Protocol # PT1200-1
Version # P936299

Protocol # PT1200-1
Version # P936299

www.clontech.com

CLONTECH Laboratories, Inc.

MTN® Blot User Manual

V. Hybridization of Oligonucleotide & cDNA Probes

For hybridizing radioactively-labeled probes follow Section A. For hybridizing nonradioactively-labeled probes follow Section B.

A. Hybridization of radioactively-labeled probes

We recommend the following probe concentrations:

- cDNA probes: 2-10 ng/ml or 1-2 $\times 10^6$ cpm/ml.
- Oligonucleotide probes: 20-50 ng/ml or 1-2 $\times 10^7$ cpm/ml.

Note: Higher probe concentrations will reduce hybridization time, but may increase background.

1. Warm ExpressHyb Solution at 68°C, and stir well to completely dissolve any precipitate. For oligonucleotide probes, equilibrate ExpressHyb at 37°C.
2. Prehybridize membranes in a minimum of 5 ml of ExpressHyb Solution, with continuous shaking for 30 min at the appropriate temperature:
 For cDNA probes: 68°C
 For oligonucleotide probes: 37°C

Note: If you are using hybridization bottles, make sure that the marked side of the membrane is flush against the side of the bottle. Bubbles between the membrane and the bottle can give the appearance of bubbles on the blot.

3. Denature radioactively labeled probes at 95-100°C for 2-5 min. Then chill quickly on ice.
4. Add radiolabeled probe to 5 ml of fresh ExpressHyb, and mix thoroughly.
5. Replace the ExpressHyb Solution with the fresh solution containing the radiolabeled probe. Remove all air bubbles from the container, and make sure ExpressHyb Solution is evenly distributed over the blot.
6. Incubate with continuous shaking for 1 hr at the appropriate temperature:
 For cDNA probes: 68°C
 For oligonucleotide probes: 37°C

7. Rinse the blot in Wash Solution 1 several times at room temperature. Wash for 30-40 min with continuous agitation; replace the wash solution several times.
8. Wash the blot two times in Wash Solution 2 with continuous shaking for 40 min at the appropriate temperature:
 For cDNA probes: 50°C
 For oligonucleotide probes: room temperature

9. Remove the blot with forceps and shake off excess wash solution.

Note: Do not allow the membrane to even partially dry. Allowing the membrane to dry can cause high background and will make subsequent probe removal difficult.

CLONTECH Laboratories, Inc.

www.clontech.com

Protocol # PT1200-1
Version # P936299

Protocol # PT1200-1
Version # P936299

www.clontech.com

CLONTECH Laboratories, Inc.

V. Hybridization Protocols continued

10. Immediately cover the blot with plastic wrap. Mount on Whatman 3 MM Chromatography paper. Wrap again with plastic wrap.
11. Expose the MTN Blot using a phosphorimaging screen. The Storm® PhosphorImager (Molecular Dynamics) is suitable for this application. Alternatively, expose to x-ray film at -70°C with two intensifying screens.
12. Strip probe from the blot by incubating the blot in sterile H₂O containing 0.5% SDS as outlined below.
 - a. Heat the sterile H₂O/0.5% SDS solution to 90-100°C.
 - b. Remove plastic wrap from blot and immediately place in the heated solution. Make sure that exposure to air is minimal.
 - c. Incubate for 10 min, shaking frequently.
 - d. Allow the H₂O to cool for 10 min before removing the blot.
 - e. Remove the blot and air-dry until it is dry enough to be slipped into a plastic bag. The membrane can be stored at -20°C until needed.

Name: Jagathpal Sheth

Date: 10/16/19

Experiment: Northern Blot analysis

16

Probe was purified in 0.5 ml of
clustering buffer, denatured at 95°C & chill mic.

Added to 7 ml of the
alkaline lysis the membrane (from

① A. Mandel, - RTN - Blot from Clontech
- used once, stopped dried) was
incubated with 7 ml of probe soln
at 68°C in a plastic bag (sealed) for $1\frac{1}{2}$ hrs.
(carefully Sheth 15)

② Purified probe added to 7 ml of
exp. hyp. ~~and~~. The plastic bag was
emptied and the ~~solution~~ bag
filled with the solution, sealed carefully
and incubated at 68°C for $1\frac{1}{2}$ hrs.

③ Discard the exp. hyp. soln and wash
place it on a dish and wash
several times with wash bf 1.
(2x SSC, 0.05% SDS) and incubate
with the same for 40 min. Replace
the wash soln 3 times (temp. RT)

④ Replace with wash soln. 2.
(0.1x SSC, 0.1% SDS, 50°C)

40 min - 3 changes.

⑤ Take the Blot in little amount of

Name: Jagathpa . Shetty

Date: 12/16, 19

Experiment:

Northern Blot analysis of CS8


17

wash buffer. Place it on a platform made out of a Whatman paper and a saran wrap.

Place another saran wrap on the top of the ~~set~~ blot immediately (do not allow it to be wet). Place the marker spots on the edges.

Place the blot inside the cassette.

with intensifying screen. expose for 18 hrs.

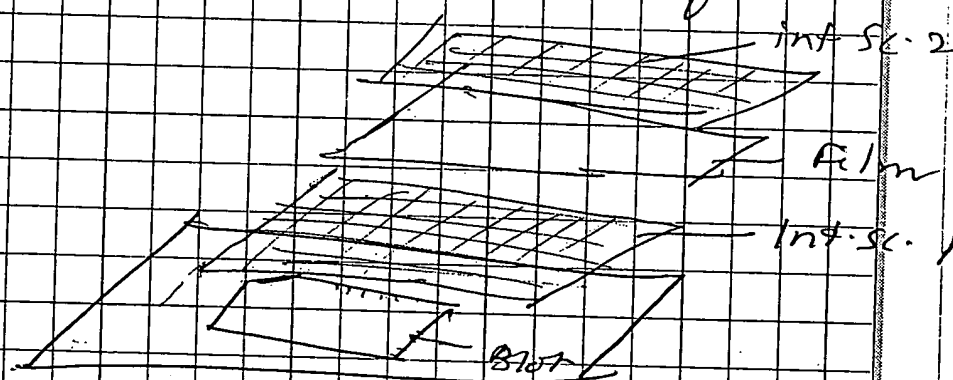

Important: Place card under blue copy. EXHIBIT 67

Name:

Jagathpai, Shekhi

Date:


12/18/19

Experiment:

Northern blot analysis - c58

18

MULTIPLE TISSUE NORTHERN
CLONTECH
developed on 12-17-99.

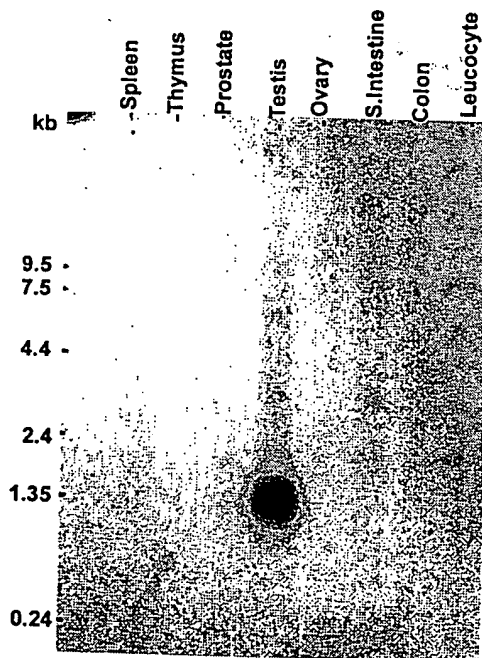

c58 is expressed only in
testis!

Important: Place card under blue copy.

EXHIBIT 68

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.